

Optilux 501. The cure-all solution.



HALOGEN CURING LIGHT



Powerful, ultra-reliable and versatile.

- **Compatibility.** Cures all light-curable dental materials on the market.
- **Four curing modes.** Provides the ultimate in versatility.
- **Built-in radiometer.** Accurately measures light. Can be used to test any halogen curing light.
- **Long-lasting bulb.** The 80w OptiBulb lasts 80 hours – output does not degrade over time.
- **Interchangeable light guides.** Accepts all Demetron light guides.



Explore the science of Kerr at www.KerrDental.com • 800.KERR.123




Optilux 501. When it comes to quality lights, we've got the cure.

At Demetron, our products – and our reputation – are built on a powerful tradition of highest quality curing lights in contemporary dentistry. Optilux 501 gives you the strongest, most versatile high-output halogen curing light available. *Ever.* Designed to be everything a clinician needs in a high-output light, the Optilux 501 truly is a cure-all solution.

The 501 has a mode for virtually any curing assignment. Continuous. Ramp. Boost. Bleach – the choice is yours. You can also choose 10-, 20-, 30-, or 40-second increments manually. The minimum output is 850 mW/cm², and in high-output modes, the output will be in excess of 1,000 mW/cm². With the 501, you're guaranteed that it will cure all light-curable dental materials on the market. Competitive high-output halogen lights can't always make that promise.

- **Radiometer.** Four-digit readout assures accurate output.
- **Patented filtering system.** Prevents high heat generation.
- **Choice of curing modes.** Multiple selections for doctor preference.
- **Programmable audible tones.** Hi-lo tone volume control is user-friendly to appeal to all preferences.
- **Bulb life indicator.** Knowledge of elapsed hours.

Optilux 501 Curing Light at 2 mm

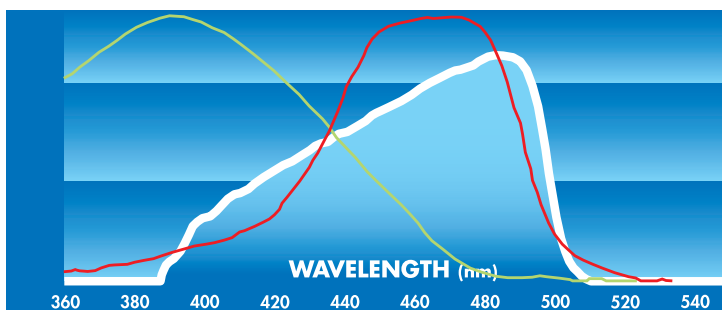

● PPD INITIATOR ● CQ INITIATOR ○ OPTILUX 501

Four Curing Modes Using the 8 mm Turbo+ Light Guide

- R Ramp.** 20-second mode with an initial 10-second ramp cure, from 100 mW/cm², and a 10-second phase in excess of 1,000 mW/cm².
- B Boost.** 10-second high-output in excess of 1,000 mW/cm².
- BLEACH Bleach.** 30-second high-output in excess of 1,000 mW/cm².
- C Continuous.** The standard curing mode, displaying time in 1-second increments. Output will be in excess of 850 mW/cm².



2006
REALITY
Five Star Award



REALITY'S Product of the Year for 2001, 2002 and 5-Star Award for 2003, 2004, 2005 and 2006

Order Information

- 910070 **Optilux 501 Curing Light**
Complete with 8 mm Turbo+ Light Guide, 11 mm Light Guide, Protective Light Shield and Hardness Disk
- 952213 **8 mm Turbo+ Light Guide**
- 20898 **11 mm Light Guide**
- 21237 **Replacement 80w OptiBulb**
- 20816 **Protective Light Shield**
(fits all manufacturers' light guides)
- 20399 **Hardness Disks** (pkg of 6 disposable)



SYBRON DENTAL SPECIALTIES

Kerr Corporation 1717 West Collins Avenue, Orange, CA 92867 800.KERR.123 714.516.7400 Fax: 714.516.7633 KerrDental.com

Lit. No. 33583